


NAVITOR™ VISION* VALVE

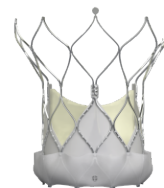

REFERENCE DIMENSIONS


Model/Reorder Number	Annulus Diameter	Annulus Area	Ascending Aorta Diameter	Annulus Perimeter
NVRO-23	19 - 21 mm	277 - 346 mm ²	26 - 36 mm	60 - 66 mm
NVRO-25	21 - 23 mm	338 - 415 mm ²	28 - 38 mm	66 - 73 mm
NVRO-27	23 - 25 mm	405 - 491 mm ²	30 - 40 mm	72 - 79 mm
NVRO-29	25 - 27 mm	479 - 573 mm ²	32 - 42 mm	79 - 85 mm
NVRO-35	27 - 30 mm	559 - 707 mm ²	27 - 44 mm	85 - 95 mm

Model/Reorder Number	Minimum Vessel Diameter	Implant Target Depth Below Annulus	Valve Inflow Diameter** A	Valve Outflow Diameter** B
NVRO-23	≥ 5.0 mm	3 mm	23 mm	23 mm
NVRO-25	≥ 5.0 mm	3 mm	25 mm	25 mm
NVRO-27	≥ 5.5 mm	3 mm	27 mm	27 mm
NVRO-29	≥ 5.5 mm	3 mm	29 mm	29 mm
NVRO-35	≥ 5.5 mm	3 mm	35 mm	35 mm

Model/Reorder Number	Aortic Stent Diameter** C	Commissure Height** D	Half Cell Height** E	Active Sealing Cuff Height**† F
NVRO-23	41 mm	21 mm	7 mm	9 mm
NVRO-25	43 mm	23 mm	7 mm	9 mm

NVRO-27	44 mm	24 mm	8 mm	10 mm
NVRO-29	46 mm	25 mm	8 mm	10 mm
NVRO-35	48 mm	27 mm	9 mm	11 mm

Model/Reorder Number	Stent Height** G	Marker Height H
NVRO-23	47 mm	3 mm
NVRO-25	48 mm	3 mm
NVRO-27	48 mm	3 mm
NVRO-29	48 mm	3 mm
NVRO-35	47 mm	3 mm

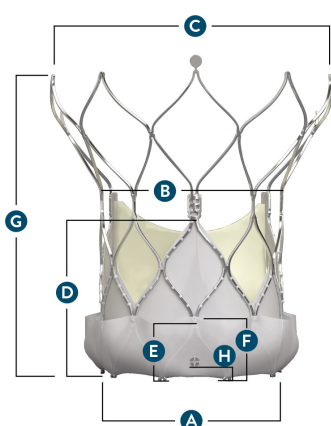

MR SAFETY INFORMATION

Non-clinical testing has demonstrated that the Navitor™/Navitor Titan™ valve (35 mm) is MR Conditional.

A patient with this device implanted can be safely scanned in an MR system meeting the following conditions:

- Static magnetic field of 1.5 Tesla (1.5T) or 3.0-Tesla (3.0T)
- Maximum spatial gradient field of 20 T/m (2000 G/cm)
- Maximum MR system reported, whole-body averaged specific absorption rate (SAR) of <2.0 W/kg (normal operating mode)

Under the scan conditions defined above, the Navitor™ valve configurations and the Navitor Titan™ configurations are expected to produce a maximum temperature rise of 5°C after 15 minutes of continuous scanning.

In non-clinical testing, the image artifact caused by the device extends radially up to 0.5 cm from the device when imaged with a gradient echo pulse sequence in a 1.5T MR system and 0.9 cm from the device when imaged with a gradient echo pulse sequence in a 3.0T MR system.

LATEX FREE

Navitor™ TAVI System with Vision Technology products and packaging do not contain latex materials based on the materials of construction. During production at some stations, manufacturing assemblers may wear non-powdered latex gloves for employee protection and prevention of contamination. Structural Heart finished goods are tested by an independent laboratory and in accordance with ISO 10993-10: 2010 Biological Evaluation of Medical Devices – Tests for Irritation and Delayed-type Sensitivities. Compliance with this standard entails animal testing under extraction conditions that exaggerate conditions of clinical use. As this testing accounts for all specified materials, environments of exposure and inherent processing, the Structural Heart division of Abbott remains confident that no sensitization or irritation responses typical of latex proteins will occur due to constituent materials, processing, and packaging used in the manufacturing of the Navitor TAVI System with Vision Technology.

NICKEL INFORMATION

The valve stent is made from nitinol, a nickel-titanium alloy that has self-expanding properties and is radiopaque. The NaviSeal™ cuff consists of an inner and outer portion which is made from polyethylene. The NaviSeal cuff provides the sealing area for implantation and provides paravalvular leak reduction by filling any potential voids between the native valve and the stent. The leaflets are made from bovine pericardium and are sutured together into a tri-leaflet configuration on the stent frame.¹

The leaflet pericardial tissue is preserved and crosslinked in glutaraldehyde. Glutaraldehyde, formaldehyde and ethanol are used in the valve sterilization process.¹

The valve leaflets are processed using Linx™ anticalcification technology. The valve is supplied sterile and non-pyrogenic in a jar containing formaldehyde storage solution. The delivery system and loading system are sterilized with ethylene oxide gas.¹

The content provided is for information only. The administration and interpretation of metal sensitivity testing is considered the responsibility of the medical team that is treating the patient with a metallic implant.

¹Navitor™ TAVI System IFU.

DISCLAIMERS AND NOTES

*Labeled as Navitor and Navitor Titan with Vision Technology

**Dimensions at fully expanded and unconstrained stent

† Outer cuff portion

#package_label#

ACCESSORIES

Model/Reorder No.	Product models	Title	Description
FNAV-DS-SM	NVRO-23, NVRO-25	FlexNav™ Delivery System, Small	The FlexNav delivery system facilitates Navitor/Navitor Vision valve implantation. The delivery system may be inserted into the vessel via its integrated sheath or through a separate arterial introducer sheath. Equivalent Integrated Sheath Diameter: 14F Valve Capsule Outer Diameter: 6.0 mm Integrated Sheath Working Length: 30 cm Delivery System Working Length: 107 cm Minimum Vessel Diameter Requirement: ≥ 5.0 mm Compatible Guidewire: 0.035 in. (0.89 mm)
FNAV-DS-LG	NVRO-27, NVRO-29, NVRO-35	FlexNav™ Delivery System, Large	The FlexNav delivery system facilitates Navitor/Navitor Vision valve implantation. The delivery system may be inserted into the vessel via its integrated sheath or through a separate arterial introducer sheath. Equivalent Integrated Sheath Diameter: 15F Valve Capsule Outer Diameter: 6.3 mm Integrated Sheath Working Length: 30 cm Delivery System Working Length: 107 cm Minimum Vessel Diameter Requirement: ≥ 5.5 mm Compatible Guidewire: 0.035 in. (0.89 mm)

CAUTION: This product is intended for use by or under the direction of a physician. Prior to use, reference the Instructions for Use, inside the product carton (when available) or at [eLabeling](#) | Abbott or at medical.abbott/manuals for more detailed information on Indications, Contraindications, Warnings, Precautions and Adverse Events.

Illustrations are artist's representations only and should not be considered as engineering drawings or photographs. Photos on file at Abbott.

Abbott

3200 Lakeside Dr., Santa Clara, CA 95054 USA, Tel: 1 800 227 9902

™ Indicates a trademark of the Abbott group of companies.

www.structuralheart.abbott

Print date: 2 Oct 2026 at 00:49:17.

